



## **ASSAY PROTOCOL**

### **Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay**

**(96-well microplate format)**

Address: Mediomics, LLC  
815 Wenneker Drive  
St. Louis, Missouri 63124 USA

Ordering: 1-800-292-4808, Direction Code 2068  
1-314-971-3028

Tech Support: 1-800-292-4808, Direction Code 2688  
1-314-971-3026

Fax: 1-314-997-2422

Website: [www.mediomics.com](http://www.mediomics.com)

## TABLE OF CONTENTS

Introduction .....	3
Bridge-It® Fluorescence Assay Platform Design .....	3
Bridge-It® L-Tryptophan Fluorescence Assay Principle.....	4
Bridge-It® L-Tryptophan Fluorescence Assay Reagents and Storage Requirements .....	5
Bridge-It® L-Tryptophan Fluorescence Assay General Notes .....	6
Bridge-It® L-Tryptophan Assay Standards .....	7
Determination of L-Tryptophan Levels in Various Samples.....	8-9
Bridge-It® L-Tryptophan Fluorescence Assay: Data Analysis .....	10
Bridge-It® L-Tryptophan Fluorescence Assay: Performance	
• L-Tryptophan Standard Curves .....	12
• L-Tryptophan Uptake by Bacteria Growing in LB Broth .....	13
• Analysis of L-Tryptophan Levels in Biological Samples .....	13-15
Specificity of the Bridge-It® L-Tryptophan Fluorescence Assay.....	16
References .....	17
Disclaimer of Warranty .....	17
Limited Use Statement .....	17
Limitation of Liability .....	18
Product Ordering Information .....	19


## Introduction

*\*The Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay Kit is intended solely for laboratory research and development (R&D) purposes. This product has not received government regulatory approval for use in diagnosis or treatment of diseases in humans or animals, monitoring of food products, or any applications other than laboratory R&D and it should not be utilized for such purposes.*

L-tryptophan (tryptophan) is one of eight essential amino acids that must be obtained from the diet. Tryptophan serves as a key building block for the synthesis of proteins and as a precursor for various brain neurotransmitters. Tryptophan is the only known precursor that can be converted into the neurotransmitter serotonin by the body. Serotonin promotes feelings of well being and calm and thereby helps to counterbalance the physiological affects of brain dopamine and the nor-adrenaline circuits that encourage fear, anger, tension, aggression, obsessive-compulsive activities, over-eating (especially of carbohydrates), migraine headache, depression and sleep disturbances. Melatonin, a metabolite of serotonin, is a sleep promoting natural hormone made by the pineal gland. In addition, tryptophan is a key precursor for niacin (vitamin B3), a vitamin that is essential for normal respiration, metabolism, and the synthesis of sex hormones. Because tryptophan plays such a critical role in the proper balancing of metabolism, mood and sleep patterns, insufficient dietary availability of this essential amino acid can lead to serious adverse clinical consequences.

## Bridge-It<sup>®</sup> Fluorescence Assay Platform Design

The common property of all sequence-specific DNA binding proteins is their ability to bind with high affinity and specificity to a DNA duplex containing a unique nucleotide sequence - i.e., the DNA binding site for the protein. Mediomics' novel assay platform design relies on this common characteristic. A DNA duplex containing the sequence-specific DNA binding site for a given target protein is split into two DNA "half-site" duplexes each having a short single-stranded overhand. These single-stranded extensions are short enough so that in the absence of the target protein little spontaneous re-association occurs. When the target protein is present, however, its high affinity for the full-length DNA sequence will drive the re-association of the two half-site DNA duplexes. This re-association can be sensitively detected by incorporating appropriate fluorescence probes into each of the two DNA half-sites. The presence of the DNA binding protein is detected as an increase in the fluorescence signal. As shown schematically below, a simple variation of this basic platform design allows a DNA binding protein (yellow oval) to function as a sensitive biosensor for its specific ligand:

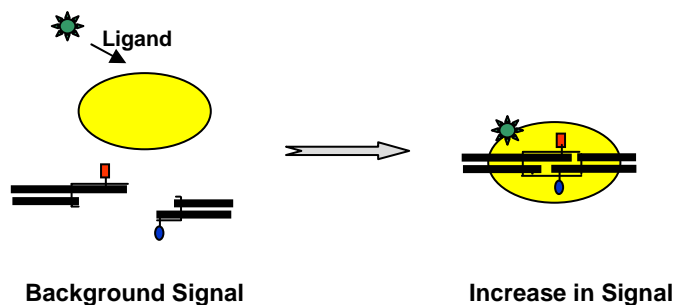

## Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay Principle

Eukaryotic cells contain an estimated 3,000 sequence-specific DNA binding proteins. These important proteins, acting either with or without a specific small molecule co-regulator (ligand), control all aspects of genomic DNA activity including gene expression, DNA replication, and DNA repair. Mediomics is applying its novel fluorescence assay platform to develop *in vitro* assays useful for rapidly and sensitively quantifying the activity of both DNA binding proteins and their small molecule ligands. L-tryptophan is a small molecule ligand that binds to and co-regulates the activity of the tryptophan repressor protein (also referred to as TrpR).

The Bridge-It<sup>®</sup> Tryptophan fluorescence assay is based on the ability of the TrpR, a bacterial DNA binding protein, to bind to its DNA binding site in a L-tryptophan dependent fashion. The two fluorochrome-labeled DNA half-fragments, each containing about one-half of the complete TrpR protein DNA binding site, are labeled one with fluorescein and the other with Oyster<sup>®</sup> 645 fluorophore. In the presence of L-tryptophan an increase in fluorescence signal can be detected as a result of the tryptophan dependent association of the two fluorophore-labeled DNA half-fragments. Tryptophan is readily detectable using the Bridge-It<sup>®</sup> tryptophan fluorescence assay in various types of test samples including bacterial growth medium, brain extract, yeast extract, as well as in human serum and urine. The linear range of the assay is 0.4  $\mu\text{M}$  - 10  $\mu\text{M}$  and the minimum tryptophan detection level is  $\sim 0.1 \mu\text{M}$ . The assay is highly specific for measuring L-tryptophan. This assay can be modified to become an assay for any enzyme reaction that uses tryptophan either as a reactant or a product.

## Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay Reagents and Storage Requirements

The quantity of reagents contained in the assay kit is adequate to perform the indicated number of individual tryptophan measurements using a black 96-well non-binding surface microplate. Store the Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay Kit in the freezer at -20°C until needed for use. **Thaw only the appropriate number of reagent tubes needed for the planned experiment.** Each tube of Tryptophan Assay Solution A contains reagent adequate for performing ten (10) tryptophan measurements using a black 96-well microplate. The tryptophan assay kit reagents will retain their activity for at least 2 months when stored frozen at -20°C. **Thaw Tryptophan Assay Solution A at 37°C for 5 minutes and keep it at room temperature until use.** Any unused reagents may be stored tightly-capped in the refrigerator at 4°C for up to one week.

Tube Cap Color	Assay Reagent	ml/tube	Storage Conditions
Green	1 mM L-tryptophan standard in water	0.25	Once thawed, store refrigerated at 4°C. <b>DO NOT RE-FREEZE.</b>
Violet	L-Tryptophan Assay Solution A	1.0	Once thawed, store refrigerated at 4°C if it is not immediately used. Stable for at least 1-week. <b>DO NOT RE-FREEZE.</b>
Blue	Buffer W*	1.0	Once thawed, store refrigerated at 4°C.

\* Buffer W may be used for diluting both the L-tryptophan standards and the test samples.

96-well black non-binding surface black microplates are recommended and individual plates may be purchased from Mediomics, LLC (Catalog # 163300)

## Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay: General Notes

The following precautionary steps are recommended in order to optimize assay performance and reproducibility:


- Pipette slowly to avoid bubbling and to allow complete recovery and transfer of assay solutions.
- If more than one tube of L-Tryptophan Assay Solution A is required for experimental analysis, mix all tubes together to ensure reagent consistency before use. Thaw only the needed number of tubes.
- After the samples have been transferred into the wells of the 96-well black microplate, gently tap the microplate on the counter top to release any trapped bubbles. Although bubbles generally break with time, large bubbles may be broken with the tip of a small gauge needle (e.g., 27 gauge) being careful to avoid loss of volume through capillary action.
- Fluorescence readings can be affected by temperature. It is recommended for consistency of assay results that all microplates used in a given study be incubated as close to the same temperature as possible (preferably ~25°C).
- To ensure consistency of results, it is recommended that a volume of 90 µl of standard or sample be transferred into the wells of a black 96-well microplate.
- High salt concentration will affect the annealing of DNA in the assay. Thus, use of high salt concentrations should be avoided in preparing the test samples and tryptophan standards to be measured using the Bridge-It<sup>®</sup> L-tryptophan fluorescence assay.
- Certain test samples (e.g., urine containing high levels of vitamins) may exhibit high background fluorescence. In order to take background fluorescence into account, it is recommended that the test sample be diluted in buffer alone (i.e., no assay solution included) and measured against a buffer blank (settings: excitation ~ 485 nm; emission: ~665 nm). Additional information concerning control of background fluorescence is presented elsewhere in this protocol (see data analysis, page 10).

## Bridge-It<sup>®</sup> L-Tryptophan Assay Standards

If more than one tube of L-Tryptophan Assay Solution A is required for experimental analysis, mix all tubes together to ensure reagent consistency before use. Thaw only the number of Assay Solution A needed for the experiment.

It is recommended that the L-tryptophan standard curve be prepared in Buffer W or under similar conditions of the test samples. If a buffer other than Assay Buffer W is used, a test to determine the effect of the alternate buffer on assay performance should be conducted prior to using it in the assay. The L-tryptophan standards should be prepared as follows:

1. Dilute the 1 mM L-tryptophan standard (green cap tube) 1:5 using Buffer W (blue cap tube) for a final concentration of 200  $\mu\text{M}$  L-tryptophan. Prepare 10 polypropylene Eppendorf tubes. Add 10  $\mu\text{l}$  of Buffer W to tubes 2-10. Add 10  $\mu\text{l}$  of the 200  $\mu\text{M}$  (2  $\mu\text{M}$ ) tryptophan stock to tubes 1 and 2. Mix the contents of tube 2 by pipetting up and down 5 times. Transfer 10  $\mu\text{l}$  from tube 2 into tube 3. Continue serially diluting the tryptophan standard through tube 9. After mixing, discard 10  $\mu\text{l}$  from tube 9. All tubes should contain a final volume of 10  $\mu\text{l}$  of tryptophan standard solution.

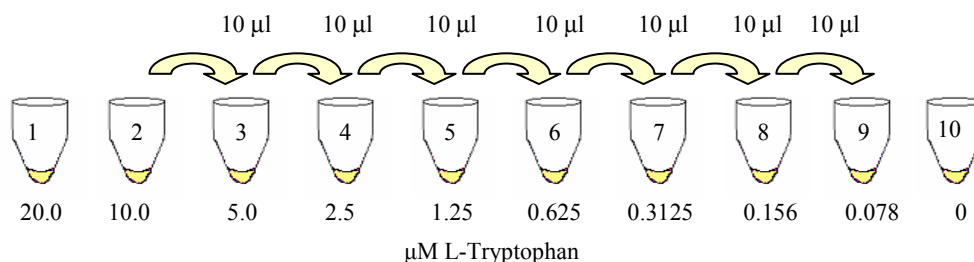

**Note:** 20  $\mu\text{M}$  is the final concentration of tube 1 when a total volume of 100  $\mu\text{l}$  is used. It is equivalent to 2 nmol/well. The L-tryptophan standards may be diluted directly in the microtiter plate. Mix by slowly pipetting solution up and down 4-5 times (avoid bubbles).

2. Add 90  $\mu\text{l}$  of L-Tryptophan Assay Solution A to each polypropylene Eppendorf tube and vortex immediately for  $\sim 1$  second at medium speed.
3. Carefully transfer 90  $\mu\text{l}$  into the appropriate well of a 96-well black microplate. The transfer of 90  $\mu\text{l}$ /well rather than 100  $\mu\text{l}$ /well when assayed directly in the microtiter plate will not affect assay sensitivity. It is recommended that the volume used for the standards and test samples be the same.
4. Cover the microplate with tinfoil to avoid evaporation and exposure to light. Incubate the microplate at room temperature ( $\sim 25^\circ\text{C}$ ) for at least 30 minutes.
5. Read fluorescence signal intensity using a fluorescence microplate reader (settings: excitation  $\sim 485$  nm; emission  $\sim 665$  nm).

## Determination of L-Tryptophan Levels in Various Samples

It is possible to use this assay to determine L-Tryptophan levels in various types of test samples including biological fluids and extracts of tissues and cells. Liquid samples (including standards) should be diluted in Buffer W. If a different buffer is used, a test to determine the effect of the different buffer on the assay performance should be conducted before using the assay. If the sample has an innate fluorescence, it is recommended a sample blank be measured at excitation ~485 nm; and emission ~665 nm (see page 10).

If more than one tube of L-Tryptophan Assay Solution A is required for experimental analysis, mix all of the Assay Solution A tubes together to ensure reagent consistency before use. Thaw only the number of tubes needed for the experiment.

### *Direct Analysis of a Liquid Sample:*

1. Clarify the sample by centrifugation if necessary. Prepare several dilutions of the sample with Buffer W to ensure that the read-out falls within the linear range of the standard curve. Add 10  $\mu$ l of sample to a polypropylene Eppendorf tube.
2. Add 90  $\mu$ l of L-Tryptophan Assay Solution A to each Eppendorf tube and vortex immediately for ~1 second at medium speed.
3. Carefully transfer 90  $\mu$ l from each Eppendorf tube into the appropriate well of a 96-well black microplate.
4. Cover the microplate with tinfoil to avoid evaporation and exposure to light. Incubate the microplate at room temperature (~ 25°C) for at least 30 minutes.
5. Read fluorescence signal intensity using a fluorescence microplate reader (settings: excitation ~ 485 nm; emission ~ 665 nm).
6. **NOTE:** It is possible to use the assay to measure up to 20% v/v of liquid sample to be added directly to the Assay Solution A with only a small affect on assay performance. To do this, add 20  $\mu$ l of the liquid test sample to 80  $\mu$ l Tryptophan Assay Solution A. The standard curve including the blank should be prepared as previously described, but an additional 10  $\mu$ l of Buffer W should be added to the 10  $\mu$ l of standard in order to bring it to the same 20  $\mu$ l volume used for the test samples. Add 80  $\mu$ l of Tryptophan Assay Solution A to each well and follow steps 4 and 5 above.

## Increase Concentration of L-Tryptophan Contained in Test Samples

If the L-tryptophan level contained in the sample extracts is found to be less than the lower limit of detection for the assay, the L-tryptophan in the sample extracts may be concentrated using lyophilization. If the protein content present in such concentrated samples is found to be problematic, ethanol extraction of the sample extract followed by centrifugation and lyophilization of the supernatant (i.e., ethanol containing L-tryptophan) is recommended. High salt concentrations affect DNA annealing and should be taken into consideration. It is not recommended that the sample extracts be lyophilized unless they are contained in a volatile solvent such as ethanol. Samples containing very low levels of L-tryptophan may be concentrated for use in the assay as follows:

1. Dilute the test samples with cold absolute (100%) ethanol to achieve a final concentration of 70% ethanol. Tightly cap all tubes, vortex and incubate on ice for 15 minutes.
2. Following the 15 minute incubation on ice, centrifuge the sample at  $\sim 12,000g$  for 10 minutes at  $4^{\circ}C$ .
3. Carefully transfer the supernatants into a fresh labeled polypropylene Eppendorf tube and lyophilize the supernatant samples (e.g., speed vac dry). Lyophilized samples that are not to be assayed for L-tryptophan content immediately should be stored at  $-20^{\circ}C$  until needed for analysis.
4. It is recommended that tubes containing the 10 assay standards (10  $\mu l$  per tube) be prepared as described on page 7, lyophilized and stored at  $-20^{\circ}C$  for assay with the lyophilized test samples.
5. Add 100  $\mu l$  of Tryptophan Assay Solution A to each of the lyophilized standards and test samples. Vortex for  $\sim 1$  second at medium speed.
6. Transfer 90  $\mu l$  into the appropriate well of a 96-well black microplate.
7. Cover the microplate with tinfoil to minimize evaporation and exposure to light. Incubate the covered microplate at room temperature ( $\sim 25^{\circ}C$ ) for at least 30 minutes.
8. Read fluorescence signal intensity using a fluorescence microplate reader (settings: excitation  $\sim 485$  nm; emission  $\sim 665$  nm).

If the tryptophan content is too high for measurement using direct dissolution in Tryptophan Assay Solution A after ethanol extraction and lyophilization, the samples may be diluted in Buffer W or in water (depending on salt concentration) for direct measurement as described above.

## Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay: Data Analysis

There are several ways of analyzing the fluorescent signal intensity readouts obtained using this assay.

1. Raw fluorescence intensity may be used as a signal proportional to L-tryptophan concentration (the x-axis is common log; the y-axis is linear).
2. Relative fluorescence (RF) change may be calculated from the raw fluorescence values. RF values are highly reproducible for the same L-tryptophan concentrations and do not depend significantly on the particular microplate reader instrument that is used to read the fluorescence signals.

RF = Relative Fluorescence  
 $F_0$  = fluorescent intensity of the 0 L-tryptophan standard (see page 7)  
 $F$  = fluorescence of L-tryptophan standard or test sample  
 $F_b$  = fluorescence of buffer blank (Buffer W only)  
 $F_{bg}$  = native sample fluorescence

- Relative fluorescence change (RF) may be calculated from the raw fluorescence values using the formula  $RF = (F-F_0)/F_0$ , or alternatively:
- Buffer and/or standard blanks may be used for the determination of a Buffer or Standard Adjusted Relative Fluorescence (Buffer Adjusted  $F-F_0/F_0$ ).
- 90  $\mu$ l of Buffer W ( $F_b$ ) is added per each of 2 wells per microtitre plate in the experiment. Read for fluorescence at the same time and conditions as the standard curves and samples. The average of the two buffer blanks is subtracted from the fluorescence readings ( $F$ ) of both the standards and samples giving a “buffer adjusted”  $F-F_0/F_0$ ” or:

$$(F-F_b) - (F_0-F_b) / (F_0-F_b).$$

- Background fluorescence present in test samples ( $F_{bg}$ ) may be accounted for by determining the fluorescence signal of the test sample versus that of Buffer W blank ( $F_b$ ) (i.e., the fluorescence of test sample in the absence of labeled probes).
- The sample (or sample dilution) should be handled in the same manner as the sample that is actually being assayed. For example, add 10  $\mu$ l of sample to 90  $\mu$ l of Buffer W ( $F_{bg}$ ), mix, and transfer 90  $\mu$ l of the solution into the well of a 96-well microplate. Read the fluorescence of the samples compared to the fluorescence readings of the buffer blank ( $F_b$ ) using a microplate reader (settings: excitation  $\sim$  485 nm; emission  $\sim$  665 nm). If  $F_{bg}$  is greater than  $F_b$ ,

subtract the F<sub>bg</sub> to eliminate the background sample fluorescence. The Relative Fluorescence (RF) then becomes


$$RF = [(F - F_{bg}) - (F_o - F_b)] / (F_o - F_b).$$

3. Calculation of the final concentration of L-tryptophan in the test sample should incorporate the dilution factor. For example, when 10 µl of the sample is used in the assay and the concentration is found to be 5 µM, the concentration of test sample is 50 µM.
4. A “Sigma Plot” software program was used to analyze and graph the results obtained from the Bridge-It<sup>®</sup> Tryptophan Fluorescence Assay. Data were converted to RF values and plotted as a sigmoidal, 4-parameter graph (x-axis as common log, y-axis as linear).


## Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay Performance

### L-Tryptophan Standard Curves

The following L-tryptophan standard curves were prepared using the Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay. Figure 1 shows the standard curve prepared in Buffer W with the fluorescence signal read-out performed after incubation for 0.5, 2, 3.5, and 5 hours after the addition of Tryptophan Assay Solution A using the 96-well format. The standard errors were from duplicate measurements. Figure 2 shows a portion (0.07-1.25  $\mu\text{M}$  L-tryptophan) of the Figure 1 standard curve read at 30 minutes and plotted to show the lower range of L-tryptophan detection sensitivity for the assay. See data analysis on page 10 for additional information.

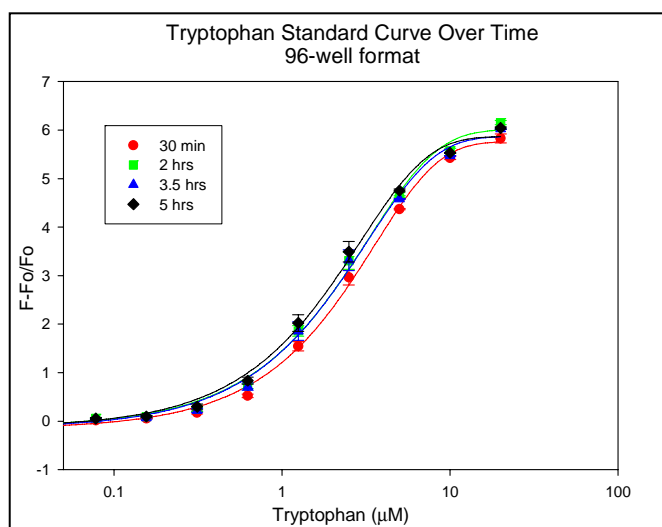

**Figure 1**

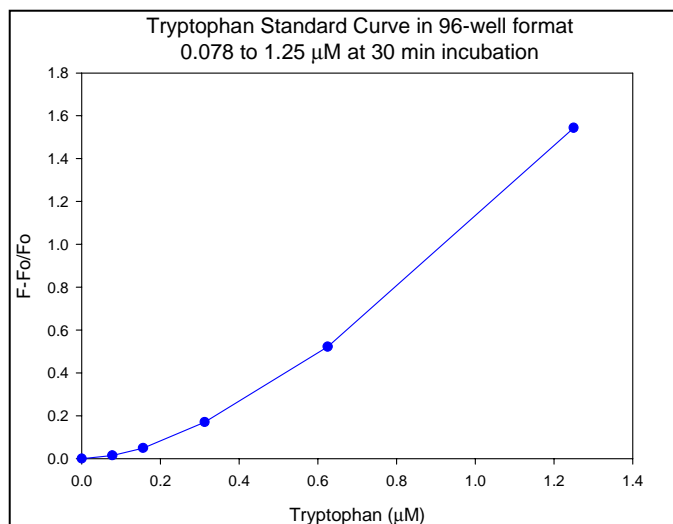


**Figure 2**

## L-Tryptophan Uptake By Bacteria Growing in LB Broth

*E. coli* strain, DH5 $\alpha$ , was seeded into LB Broth and grown overnight at 37°C. The following day, the culture was diluted to OD<sub>600nm</sub> = 0.1, divided into two flasks, and incubated at 37°C with shaking. Samples were taken each hour for a period of up to 7 hours. Bacterial growth in the duplicate cultures was measured at OD<sub>600nm</sub> at each time point.

Hourly samples were also taken and centrifuged at 12,000 rpm at 4°C for 10 minutes. The supernatants were aliquoted and immediately frozen. The frozen samples were subsequently thawed, diluted in Assay Buffer W, and analyzed for L-tryptophan content using both the Bridge-It<sup>®</sup> L-tryptophan fluorescence assay method and the AccQ-Tag HPLC procedure. Using the Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay, the LB Broth alone was found to contain 331.7  $\mu$ M L-tryptophan at the beginning of the culture. The initial concentration of tryptophan is considered 100%. The results of these studies are presented in Figure 3.

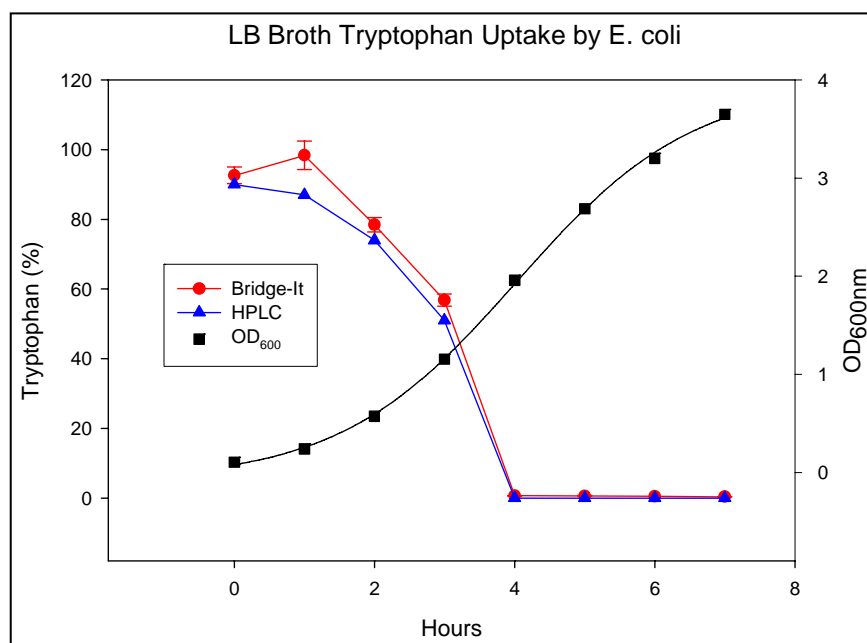

**Figure 3**

## Analysis of L-Tryptophan Levels in Biological Samples

### Rat Brain Extract:

Rat brain was extracted using the following procedure. Two volumes per gram (v/g) of brain of 50 mM  $\text{NH}_4\text{HCO}_3$  were used to homogenize the rat brain tissue. After centrifugation at 12,000g for 30 minutes at 4°C, the supernatant brain extract was transferred into another centrifuge tube and diluted with cold absolute (100%) ethanol to a final concentration of 70%. The tube was tightly capped and placed into an ice bath. After 30 minutes of incubation on ice, the ethanol-treated extract was clarified by centrifugation. The supernatant solution was aliquoted into 1 ml samples in Eppendorf tubes, lyophilized (speed vac dried) and stored frozen at -20°C until needed for analysis. An aliquot was later resuspended in 80  $\mu\text{l}$  of cold Buffer W and centrifuged at 12,000 rpm for 15 min at 4°C. Results of the analysis of the rat brain extract for the presence of L-tryptophan are shown on Table 1. The L-tryptophan level in the rat brain extract that was prepared and analyzed in our laboratory was determined to be 7.2 nmol L-tryptophan per gram of wet weight rat brain tissue.

### Yeast Cell Extract:

Yeast strain BY4743 was grown in Yeast-Bactopeptone-Dextrose (YBD) broth overnight at 30°C with gentle agitation. The following day, the culture was diluted to a density of approximately 0.02  $\text{OD}_{600\text{nm}}$  and again cultured overnight at 30°C with gentle agitation. The yeast cells were harvested at a density of 8.0  $\text{OD}_{600\text{nm}}$  by centrifugation at 4°C. The yeast cells were washed 2-times with cold PBS and the cells pellet was resuspended in  $\sim 1/30^{\text{th}}$  of the original volume of the culture using cold double distilled water. An equal volume of washed glass beads (0.5 mm) was added to the suspended cell pellet and vortexed at top speed for one minute followed by incubation on ice for one minute. This cycle of vortexing and incubating on ice was repeated 6-times. The sample was centrifuged and the supernatant yeast cell extract was diluted using absolute (100%) ethanol to a final concentration of 70% ethanol and incubated on ice for 10 minutes. A final centrifugation removed the resulting precipitate. The clarified supernatant containing the ethanol extracted L-tryptophan was divided into 1 ml aliquots, dried in a speed vac, and immediately frozen for storage in a freezer at -20°C. The following day, a frozen aliquot was thawed, reconstituted using 30  $\mu\text{l}$  cold autoclaved double distilled water, diluted 1:8 in Buffer W and centrifuged at 4°C to remove any precipitation prior to performing the tryptophan analysis using the Bridge-It<sup>®</sup> L-tryptophan fluorescence assay. The results of this analysis are presented in Table 1. The tryptophan level in the yeast cells was determined to be 83.84 nmole L-tryptophan in  $1 \times 10^9$  yeast cells.

### Human Urine:

Samples of urine from two normal adult human subjects were centrifuged and the clarified urine (supernatant) was aliquoted (1 ml/tube) into clean labeled Eppendorf tubes and immediately frozen for storage in the freezer at -20°C until needed for testing. Subsequently the frozen urine samples were thawed, diluted in cold Assay Buffer W, and centrifuged again to remove any precipitate. The clarified supernatant urine was analyzed for L-tryptophan content using the Bridge-It<sup>®</sup> L-tryptophan fluorescence assay. Urine dilutions of 1:4 or less were found to have a significant level of background fluorescence. This background fluorescence was taken into account as part of the assay data analysis. The L-tryptophan

levels found in two normal human urine specimens using the Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay are presented in Table 1.

### Human Blood:

Commercially-available, normal, pooled human serum and pooled plasma and also fresh human serum and plasma from a single normal donor were simultaneously analyzed for L-tryptophan using the Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay. Following dilution in cold Buffer W, the serum and plasma samples were centrifuged to remove any precipitate prior to L-tryptophan analysis. The L-tryptophan levels in fasting normal human plasma and serum samples have been reported in the literature to be in the range of 24-79  $\mu\text{M}$  for children (1-16 years old) and 20-95  $\mu\text{M}$  for adults<sup>5</sup>. As shown in Table 1, the levels of L-tryptophan obtained for the samples of normal human serum and plasma using the Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay fell well within the reported normal range for L-tryptophan in normal human blood samples.

**Table 1**


L-Tryptophan Levels Found In Various Types of Test Samples.

Sample	Bridge-It <sup>®</sup> L-Tryptophan ( $\mu\text{M}$ )
Rat Brain Extract	13.6 $\pm$ 0.6
Yeast Cell Extract	310.5 $\pm$ 6.6
Human Urine 1	27.6 $\pm$ 5.3
Human Urine 2	21.8 $\pm$ 7.0
Pooled Human Serum	43.7 $\pm$ 1.9
Pooled Human Plasma	35.9 $\pm$ 0.7
Fresh Human Serum	27.5 $\pm$ 3.2
Fresh Human Plasma	22.5 $\pm$ 3.6

Note: The data represents the average value  $\pm$  standard deviation (n = 4)

## Specificity of the Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay

To assess the specificity of the Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay, all 20 amino acids were tested up to 100  $\mu\text{M}$  concentration. Nineteen (19) amino acids (excluding L-tryptophan), were shown to have no cross-reactivity with L-tryptophan when they were measured using the Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay. Further, D-tryptophan (minimum of 98% pure as measured by thin layer chromatography), serotonin, and 5'-HTP (a precursor of L-tryptophan) were also analyzed up to 20  $\mu\text{M}$  concentration using the Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay and no cross-reactivity observed as demonstrated by the results shown below in Figure 4.

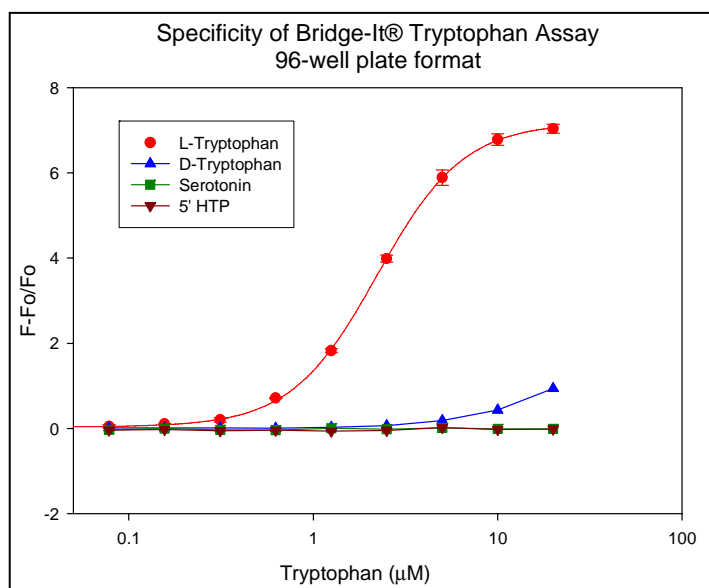

**Figure 4**

## References

1. Heyduk, T. & Heyduk E., Molecular beacons for detecting DNA binding proteins. *Nature Biotechnology* **20**, 171-176 (2002).
2. Heyduk, E., Fei, Y., & Heyduk, T., Homogeneous Fluorescence Assay for cyclic AMP, *Combinatorial Chemistry and High-Throughput Screening*, **6**, 183-194 (2003).
3. US Patent No. 6,544,746, Rapid Proximity-Based Assay for the Detection and Quantification of DNA-Binding Proteins.
4. Heyduk T, Knoll E, & Heyduk E, Molecular beacons for detecting DNA binding proteins: mechanism of action. *Anal. Biochem.* **316**, 1-10 (2003).
5. Clinical Guide to Laboratory Tests, Edited by Norbert W. Tietz, 3<sup>rd</sup> Edition, W.B. Sanders Company (1995).

## Disclaimer of Warranty

Despite Mediomics best efforts to provide complete and accurate information to its customers, it is not possible to completely insure that all information provided is current and accurate. The information contained in this document could potentially include technical inaccuracies or typographical errors. Mediomics does not assume responsibility or liability for any actions taken as a result of using Mediomics' products. Mediomics assumes no responsibility for errors or omissions in the content included in this product. **THESE ARE PROVIDED "AS IS" WITHOUT ANY WARRANTY OF ANY KIND, EITHER EXPRESS OR IMPLIED, INCLUDING BUT NOT LIMITED TO IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, OR NON-INFRINGEMENT OF THIRD PARTY RIGHTS.** Mediomics does not warrant the accuracy or completeness of text, graphics, or other information contained herein. In no event shall Mediomics be liable for any special, indirect, incidental or consequential damages of any kind, or any damage whatsoever resulting from loss of use, data, or profits, whether or not advised of the possibility of damage, and on any theory of liability, arising out of or in connection with the use of these materials. All information, products and services that are referenced to in this Protocol are provided subject to all applicable terms and conditions. Mediomics provides this Protocol as a service to its customers and the public solely for informational purposes. Reproduction of any of the content contained in this document is prohibited.

## Limited Use Statement

This product is the subject of US Patent 6,544,746 which is exclusively licensed by Mediomics, LLC from Saint Louis University. The purchase of this product conveys to the buyer the non-transferable right to use the purchased amount of the product and components of the product solely for laboratory research conducted by the buyer (whether the buyer is an academic or for-profit entity). The buyer cannot sell or otherwise transfer (a) this product (b) its components or (c) materials made using this product or its components to a third party or

otherwise use this product or its components or materials made using this product or its components for Commercial Purposes. The buyer may transfer information or materials made through the use of this product to a scientific collaborator, provided that such transfer is not for any commercial purpose, and that such collaborator agrees in writing (a) to not transfer such materials to any third party, and (b) to use such transferred materials and/or information solely for research and not for Commercial Purposes. The term "Commercial Purposes" as used herein means any activity by a party for consideration and may include, but is not limited to: (1) use of the product or its components in manufacturing; (2) use of the product or its components to provide a service, information, or data; (3) use of the product or its components for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the product or its components, whether or not such product or its components are resold for use in research. Mediomics, LLC will not assert a claim against the buyer for infringement of the above patents based upon the manufacture, use, or sale of a therapeutic, clinical diagnostic, vaccine or prophylactic product developed in research by the buyer in which this product or its components was employed, provided that neither this product nor any of its components was used in the manufacture or manufacturing process quality control of such product. If the purchaser is not willing to accept the limitations of this Limited Use Statement, information on purchasing a license for use of this product for purposes other than laboratory research, contact Mediomics, LLC, Business Development, 815 Wencker Drive, St. Louis, MO 63124, Phone: (314) 971-3028, Fax: (314) 997-2422.

### **Limitation of Liability**

Mediomics shall not be liable for any damages suffered as a result of using, modifying, contributing, copying, or distributing any of the information contained in this document. In no event shall Mediomics be liable for any indirect, punitive, special, incidental or consequential damage (including loss of business, revenue, profits, use, data or other economic advantage), however it arises, whether for breach or in tort, even if Mediomics was previously advised of the possibility of such damage. User has sole responsibility for adequate protection and back-up of data and/or equipment used in connection with use of the Mediomics' website and documentation and will not make any claim against Mediomics for lost data, rerun time, inaccurate output, work delays or lost profits resulting from the use of the materials contained in this product. Moreover, the user agrees to hold Mediomics harmless from, and covenants not to sue Mediomics for, any claims based upon the handling and use of this product by the user.

## Product Ordering Information

### Bridge-It<sup>®</sup> L-Tryptophan Fluorescence Assay Kit (96-well microplate format)

Product Description	Measurements/ Kit	Catalog No.	Price <sup>3</sup>
Bridge-It <sup>®</sup> L-Tryptophan Fluorescence Assay Kit <sup>1,2</sup>	50	1-1-1001A	\$164.00
Bridge-It <sup>®</sup> L-Tryptophan Fluorescence Assay Kit <sup>1,2</sup>	96	1-1-1001B	\$315.00

### 96-well black non-binding surface microplate:

Product Description	Catalog No.	Price <sup>3</sup>
96-well black non-binding surface microplate	163300	\$6.25

<sup>1</sup> Bridge-It<sup>®</sup> is a registered trademark of Mediomics, LLC, St. Louis, Missouri, USA.

<sup>2</sup> Oyster<sup>®</sup> is a registered trademark of Denovo Biolabels, GmbH, Munster, Germany.

<sup>3</sup> All prices are denominated in U.S. dollars. Shipping and handling costs will be applied. Prices shown may be changed without notice.

- **To Order:**

E-mail - [orders@mediomics.com](mailto:orders@mediomics.com)

Telephone -

Toll free: 1-800-292-4808, Direction Code 2068;

Direct dial: 1-314-971-3028

Fax: 1-314-997-2422.

- **For Customer Service:** E-mail - [customerservice@mediomics.com](mailto:customerservice@mediomics.com)

Telephone -

Toll-Free: 1-800-292-4808, Direction Code 2068

Direct Dial: 1-314-971-3028

- **For Technical Support:** E-mail- [techsupport@mediomics.com](mailto:techsupport@mediomics.com);

Telephone -

Toll-Free: 1-800-292-4808, Direction Code 2688

Direct Dial: 1-314-971-3026